




**verilux**<sup>®</sup>  
The Healthy Lighting Company<sup>®</sup> Est. 1956

**Eliminate up to 99%  
of germs, viruses and  
bacteria from your  
toothbrush.**

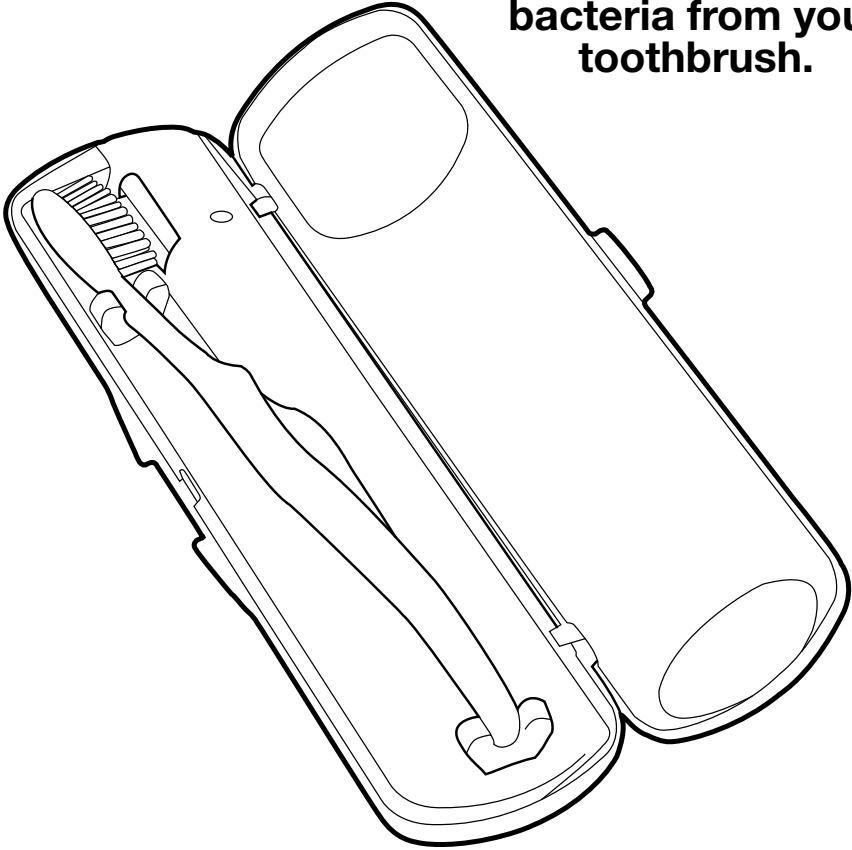

**CleanWave** <sup>®</sup>  
Portable Toothbrush Sanitizer

**VH17**



Dear Customer,

Thank you for purchasing the Verilux CleanWave® Portable Toothbrush Sanitizer. You have received an innovative and quality manufactured product, backed by a one-year limited warranty.

Many other healthy lighting products are available through our catalog and online. Visit us on the web at [www.verilux.com](http://www.verilux.com) to request a free catalog and learn more about all our quality Verilux products, or call our toll-free number 1-800-454-4408.

As a Verilux customer, your satisfaction means everything to us. We look forward to serving you now and in the future.

Have a Bright Day!

Nicholas Harmon  
President, Verilux, Inc.

## Precaution

- Do not mix old and new batteries.

## Warnings

- Do not immerse in water or liquid of any kind.
- UV-C light is emitted from the light. Avoid looking directly into the illuminated light when the sanitizer is turned on.
- The bulb emits UV-C light. If the bulb remains lit when the lid is open, do not look directly into the bulb. Remove the batteries and contact customer service.
- Dispose of device per local regulations.
- Keep out of reach of children under 6 years old.

## Safeguard

- Unit is equipped with a safety mechanism that should not allow it to go on when the lid is open.

## Features

- Auto start with built-in timer and safety switch.
- UV-C light eliminates up to 99% of viruses, germs and bacteria.
- Sanitizes manual brushes and electric toothbrush heads.
- Sturdy case keeps your toothbrush clean while traveling or in the bathroom.

## Operation

Operating your Verilux® CleanWave® Portable Toothbrush Sanitizer:

1. Turn Sanitizer over and remove the lid on the battery compartment.
2. Insert 2-AAA alkaline batteries (not included), making sure they are inserted in the proper direction.
3. Replace the battery compartment cover.
4. Open the lid of the Sanitizer and insert your toothbrush with the neck in the clip (make sure the bristles are facing the light at the back of the compartment).
5. Close the lid, and you will notice the light has illuminated (the lid will protect your eyes from the UV-C light).
6. The built-in timer will provide a complete 6-minute sanitizing cycle before automatically turning the Sanitizer off.
7. Although the unit has air circulation holes, it is recommended that you shake off all excess water prior to placing your toothbrush in the Sanitizer.
8. There is no need to wait after the sanitizing cycle is complete to use your toothbrush again.

## Technical Specifications

### Environment

Operating Temperature: 15°C to 35°C (59°F to 95°F)

Storage Temperature: 0°C to 40°C (32°F to 104°F)

Relative Humidity: 5% to 95% non-condensing

### UV-C Bulb

0.6 Watt UV-C bulb; non-replaceable

### Battery

2-AAA alkaline batteries (not included)



## One Year Limited Warranty

**ATTENTION! ONCE OPENED, PLEASE DO NOT RETURN THIS PRODUCT TO THE STORE WHERE IT WAS PURCHASED FOR REPAIR OR REPLACEMENT!**

Many questions may be answered by visiting [www.verilux.com](http://www.verilux.com), or you may call our Customer Service Department at 800-786-6850 during normal business hours.

This limited warranty is provided by: Verilux, Inc., 340 Mad River Park, Waitsfield, VT 05673

Verilux warrants this product to be free from defects in material and workmanship for a period of one year from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, Inc. will, at its option, repair or replace defective parts of this product, at no charge to the customer, subject to the following limitations: This limited warranty does not include any postage, freight, handling, insurance or delivery fees. This warranty does not cover damage, defect or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. The use of shipping insurance is recommended to help protect your product.

**Return Authorization is required for all returns.** To obtain a Return Authorization, please contact the Verilux Customer Service Department at 800-786-6850.

If, during the first year of ownership, this product fails to operate properly, it should be returned as specified at [www.verilux.com/warrantyreplacement](http://www.verilux.com/warrantyreplacement) or as instructed by a Verilux customer service representative at 800-786-6850.

**Note:** *Verilux recommends using a quality surge suppressor on all electronics equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges and may be purchased at electronics stores.*

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

**Please visit our website at [www.verilux.com](http://www.verilux.com) or call our Customer Service Department at**

**1-800-786-6850**

Representatives are available Monday — Friday, 9am — 5pm EST

**verilux®**

340 Mad River Park, Waitsfield, VT 05673

340 Mad River Park, Waitsfield, VT 05673